CA15FL

Fuel Cell Breathalyzer

Owner's Manual – Multilanguage Edition

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CA15FL

Owner's Manual


Fuel Cell Technology
Digital Alcohol Detector

This device is a breath alcohol tester designed to measure the presence of alcohol in the body. Please make sure to read this manual thoroughly before use.

DO NOT DRINK AND DRIVE


1. Part Names

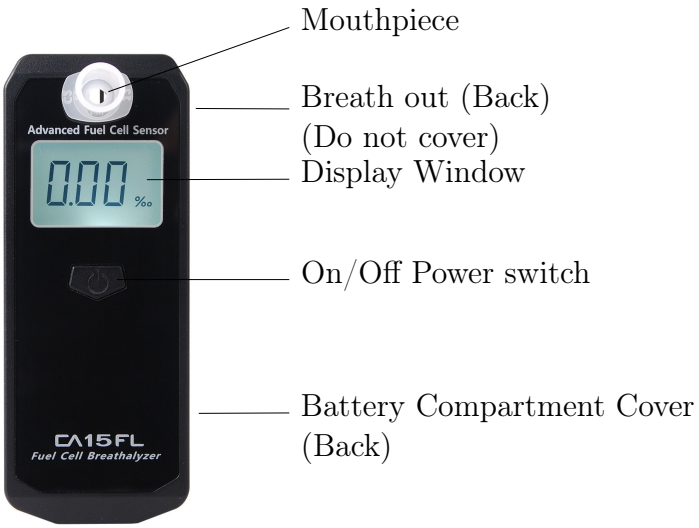- **Mouthpiece** – insert a mouthpiece here and blow into it.
- **Breath out (Back)** – exhaust on the back of the unit; do not cover it.
- **Display Window** – shows the test result and status signals.
- **On/Off Power Switch** – powers the unit on and off.
- **Battery Compartment Cover (Back)** – covers the two AAA batteries.


2. Battery Replacement

Apply a small amount of pressure on the battery cover located on the back of the unit and push down and away from the unit to release the cover. See the following diagram.

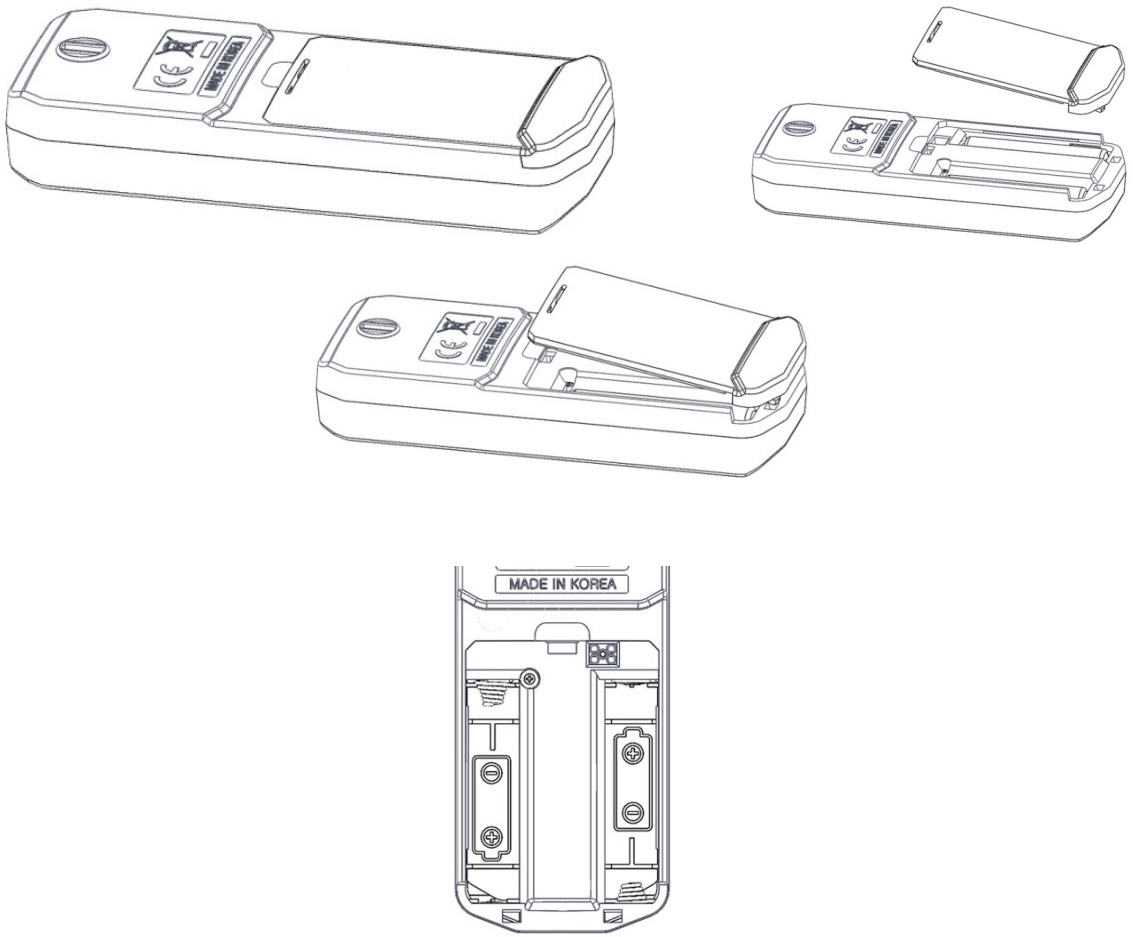∗ **N.B.:** Make sure to insert the battery with the proper position of positive (+) and negative (−) as indicated.

3. Display

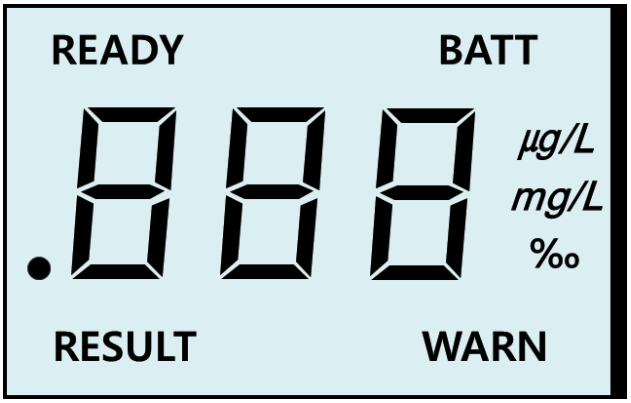- **READY** – start blowing when the READY signal is on the display.
- **BATT** – the BATT signal means low battery voltage. It is recommended to replace the 1.5V Alkaline Batteries.
- **RESULT** – the test result is displayed with the RESULT signal.
- **WARN** – if the result is over 0.50‰ (or equivalent to it in other measurement standard) the WARN signal is displayed.
- **Reading** – the test result is displayed in BrAC / Promile data, correspondent to the 3rd decimal readout at BAC.
- **Measurement unit** – the measurement unit (μg/L, mg/L, ‰) is set to the region standard.



This unit has a Dual Color LCD Backlight. See the **Error (Warning) Notification** section in the manual.

4. Operation

1. Insert a mouthpiece into the mouthpiece inlet.
2. Push the ON/OFF power switch that is located on the right side under the display window.
 - The unit will show “Full display” firstly and “on” along with a welcome sound. Then it starts a countdown which takes about 15–35 seconds. This warm-up process prepares the sensor and circuit for testing.
3. When you hear a beep sound and the “bLo” display is on, please blow steadily into the mouthpiece until you hear the “Beep” sound again. Stop blowing.
 - In this stage, if you don’t blow within 20 seconds, it automatically turns off. Try to power on again.
4. After the circling line is on the display for 3–22 seconds, the test result (mg/L BrAC or ‰ Promile) will be displayed for 20 seconds along with a beep sound.
 - If the test result is over 0.25mg/L / 0.50‰ (default setting) respectively, the LCD color will change to red and the “WARN” signal will be displayed along with an “Alarm” sound.
 - The result will be displayed for 20 seconds (light on for the first 15 seconds and light off for another 5 seconds). If you push the power switch during the last 5 seconds, the unit goes back to the countdown stage for another test.
5. Finally, the unit will show “OFF” and turn off automatically along with a sound.
6. For successive testing, try again from No.2 to No.4.
7. Space tests at least 2 minutes apart.

∗ After using, you can turn the device off by pressing the power switch for 2 seconds.

5. Error (Warning) Notification

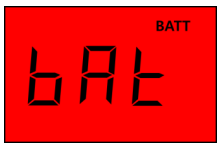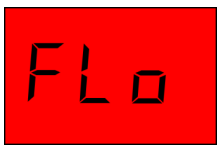Error Message	Cause & Action
	If the result is higher than 0.25mg/L BrAC (0.50‰ Promile), the "WARN" signal is displayed along with an "ALARM" sound.
	If the result is above 2.00mg/L BrAC (4.00‰ Promile), the "Hot" signal is displayed along with an alarm sound.
	When the test count is over 500, the unit shows "CAL". Even if the unit still works, re-calibration is required. Consult the seller or local distributor for the service.
	The battery voltage is too low to perform a test. The device will turn off automatically after showing this signal. Replace the 1.5V Alkaline Batteries.
	The user had a short exhale time during the test, which likely means the user did not exhale steadily and continuously for at least 3–5 seconds. Push the power button to perform another test.
	When the ambient temperature is out of the working temperature range (5 ℃ – 40 ℃), the unit shows "OUT" and then turns off. Keep and use the unit within the recommended temperature range.

∗ **N.B.:** When any of these errors occurs, the LCD Backlight turns red to give notice.

6. Precautions

1. After drinking, smoking and eating, it is recommended to test after 20 minutes.
 - It takes approximately 20 minutes for alcohol to be absorbed into the blood from the digestive organs, and residual alcohol remaining in the mouth takes this long to dissipate.
2. Avoid testing under strong wind or in a closed room with contaminated air.
3. When the "bAt" signal is on the display, it is advisable to replace the 1.5V alkaline batteries.
4. Do not blow cigarette smoke, food or liquid into the instrument, because this may damage the sensor.
5. Do not use this unit as a tool to determine whether you or anyone else should operate any motor vehicle or machinery.
6. Keep it under 0–40 ℃ (32–104 °F). This device is designed to be used in a temperature range of 5–40 ℃ (41–104 °F). Keep and use the unit within the recommended temperature range.
7. Correlation between breath alcohol content and blood alcohol concentration depends on many variables, including environmental factors (such as air quality, wind, humidity, temperature, etc.) and the health conditions of the test subject.
8. In normal use, the CA15FL requires re-calibration every 6 months or after 500 tests to obtain accurate results. Discuss with local distributors for re-calibration.
9. You cannot use the results of this unit in a court of law. The driver holds full responsibility for driving under the influence of alcohol.
10. Do not use this device in conjunction with an ignition interlock device or while sitting in a motor vehicle.
11. An individual's testing result can continue to increase or decrease over time, depending on many factors. Obtaining a testing result reading at one time does not mean the testing result reading will be the same even a few minutes later.
12. Do not let yourself or anyone else drink and drive. Always use a designated driver when you or anyone else drinks alcohol.

7. Specifications

Sensor	Advanced Fuel Cell Sensor
Accuracy	±0.005% BAC at 0.050% BAC
Detection Range	0.00 – 2.00mg/L BrAC 0.00 – 4.00‰ Promile
Warm-up Time	15 – 35 seconds
Response Time	3 – 22 seconds
Dimension	111.5mm × 45mm × 20.2mm
Weight	58.5g (81.2g battery included)
Power Supply	Two 1.5V "AAA" size Alkaline Batteries
Working Temp.	5 ℃ – 40 ℃
Storage Temp.	0 ℃ – 40 ℃
Calibration	Every 6 months or after 500 tests
Warranty	1 Year Limited Warranty

8. Warranty

The manufacturer warrants for 1 year that the product is free from defects in workmanship or material (excluding calibration) in normal service from the date of purchase.

The manufacturer's obligation under the warranty is limited to replacing, adjusting, or repairing the unit if returned along with proof of purchase.

This warranty is void if the unit has been tampered with or maliciously damaged.

Distributor
Sunnysoft s.r.o.
Kovanecka 2390/1a
190 00, Praha 9 - Liben
www.sunnysoft.cz

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